

PROTOCOL: FINGER WO ROUTINE

COIL:	Dedicated hand/wrist coil. If not present, small flex coil (4-16 CH) wrapped around.								
POSITIONING:	Fingers straight. Pad fingers dorsal and volar so they cannot move. Superman arm at side. If arm at side careful attention to avoid phase wrap. Extra effort make patient comfortable.								
CLINICAL INDICATIONS/ HISTORY:	FINGER PAIN, FRACTURE, LIGAMENT TEAR, TENDON TEAR, ARTHRITIS, INJURY, CONTRACTURE, VOLAR PLATE TEAR								
PLOTTING:	Axial FIRST off of coronal and sagittal localizers. True axial to finger Coronal = coronal to finger of interest = sagittal to finger of interest								
COVERAGE:	May vary depending on area of involvement. Approx fingertips to metacarpals								
TIPS:	Note that finger motion causes artifactual signal in tendons. Difficult to hold fingers still for entire scap. Extra padding to hold still.								

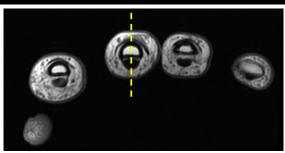

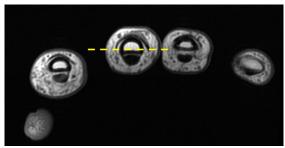

SCAN ORDER	PLANE	IMAGE CONTRAST/ WEIGHTING	MODE/ SEQ TYPE	Averages (NEX)	BW (Hz/pixel)	ETL	TR RANGE	TE RANGE	ті	FLIP ANGLE	SLICE/GAP (mm)	FOV (cm)	Resolution	Phase Axis	Send to PACS	Fat sat
1	АХ	T2FS	2D/FSE	1	150-180	8-12	3500- 6000	65-75		>130	3.0/0	8-12	256 x 256	RL	FULL SERIES	WEAK
2	АХ	PD	2D/FSE	2	150-180	8-12	3500- 6000	38-45		>130	3.0/0	8-12	512 X 320	RL	FULL SERIES	NONE
3	COR	PDFS	2D/FSE	1	150-180	8-12	3500- 6000	38-45		>130	1.5-2.0	8-12	256 X 256	SI	FULL SERIES	WEAK
4	COR	PD	2D/FSE	2	150-180	8-12	3500- 6000	38-45		>130	1.5-2.0	8-12	512 x 320-384	SI	FULL SERIES	NONE
5	COR	T1	2D/FSE	1	150-180	2-3	450-650	8-10		>130	1.5-2.0	8-12	256 x 192	SI	FULL SERIES	NONE
6	SAG	PD	2D/FSE	1	150-180	8-12	3500- 6000	38-45		>130	2.5/0	8-12	512 x 320	SI	FULL SERIES	NONE